

3


Maintenance Checklist

Operating procedure for order no.

System components to be inspected:

- ☐ ID / PRID - Station SN..... SN.....
☐ ADC QS Server Station SN
☐ ADC QS Client Station SN

Device	Test Point - Maintenance Point	OK	not OK
Complete System	Talk to the operators about possible problems. Talk to the physician in charge about the image processing parameters in the various menus.		
ADC QS – Server Station	Clean the ventilation openings of the Workstation. Clean the interior fans.		
ADC QS Server Station and Periphery	Check all external cables and cable connections for damage and loose contacts.		
ID – Station	Clean the PC (ventilation openings).		
	Check the function of the ID Station by identifying a cassette.		
ID – Tablet	Clean the surface of the ID tablet.		
	Check the cassette duct for tight mounting.		
ID – Console	Check the console for damage.		
	Check all cables and cable connections for damage and loose contacts.		
	Check that all protection clips of the tables are present.		
	Check the protection (cable ties) of the monitor.		
Monitor	Check calibration of the monitor by means of a SMPTE test pattern. If an image error was detected, talk to the radiologist to decide the recalibration of the monitor. Recalibrate the monitor or contact the responsible Application Specialist.		
	Clean the monitor.		
Complete system – Image quality	Expose a grey field (flatfield) and evaluate the image on the workstation and the Laserprinter. Criteria: homogenous field, stripes a.s.o. Recalibrate if necessary.		

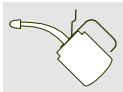

Device	Test Point - Maintenance Point	OK	not OK
Complete system - Image quality	Check plates for running askew with a test sheet. Check test sheet for wobbling of the Galvo in digitizer.		
Maintenance completion	Enter the <Maintenance date> in the service program.		
	At the end of the maintenance, hand over the maintenance checklist to the customer and inform him about the results of the maintenance.		
	If any image faults were detected, inform the responsible Application Specialist.		

Remarks:

.....

.....

.....
Date / Signature Service Technician.....
Customer